



Manual

Chromotherapy equipment with high brightness RGB LED. It incorporates a complete range of colors and different emission possibilities, with automated and extremely easy to use programs.




Material: Aluminum
 Light: 25-30 lumen
 Power: 3 AAA alkaline batteries

- 1- Power button: ON/OFF.
- 2- Function button.
- 3- Unscrew cover for changing batteries


FUNCTIONS:

- Press once to change the color of the light.
- After going through the 10 different colors, you will enter sequential mode:
 - * After magenta:
 - 1^a Sequential program press 1 colour per second.
 - 2^a progressive sequential program.

Sequence:

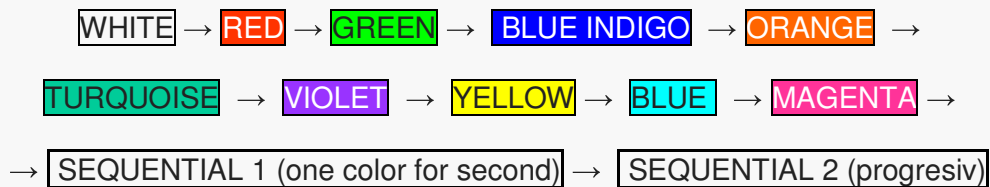

- Press for 3 seconds to activate the flashing mode (light pressed):
 Press once at 2.5 Hz. - 2 times press at 0.7 Hz. - 3 times to return to continuous broadcast mode.
- Press for 4 seconds to pause, it will activate in the last color when you press the button again. press progressive sequential program.

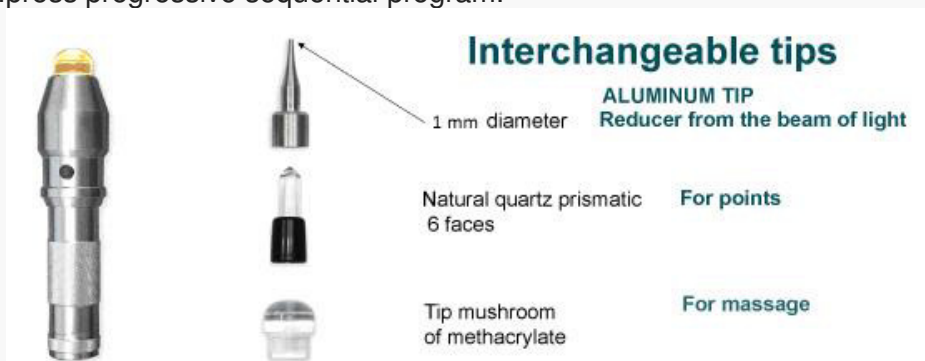

Consell de Cent, 146 08015 Barcelona T 934 120 222 F 934 121 211
www.zenlong.es zenlong@zenlong.es

